




You refine the best research works; we offer the best lighting options to assist you achieve goals.



YODN, a special lighting company, proudly announces the introduction of the LED illuminator product, designed specifically for fluorescence microscopes, YODN Hyper E301. The illuminator equips the best LED light engine, produces spectrums that can fulfill most of the research lighting demands with its wide-field spectrums design methodology. Regardless you are a researcher or a microscopist, Hyper E301 helps you conduct your works easily in an energy efficient way. YODN Hyper E301 illuminator successfully adapts the latest LED technology, remarkably produces spectrums from 350 nm to 700 nm, a true wide-field illuminator that can support different fluorescence excitation, such as DAPI, GFP/FITC, mCherry/Texas Red, Cy3, and more. Both manual and automatic digital programming control are seamlessly integrated in Hyper E301, offer the users easy, yet flexible operation options. Better yet, the LED technology long life characteristic makes Hyper E301 an almost maintenance free equipment, no more light bulb changes. YODN Hyper E301 is the best choice for fluorescence microscope applications.

HYPER E301

Wide-Field Excitation Illuminator for Fluorescence Microscope

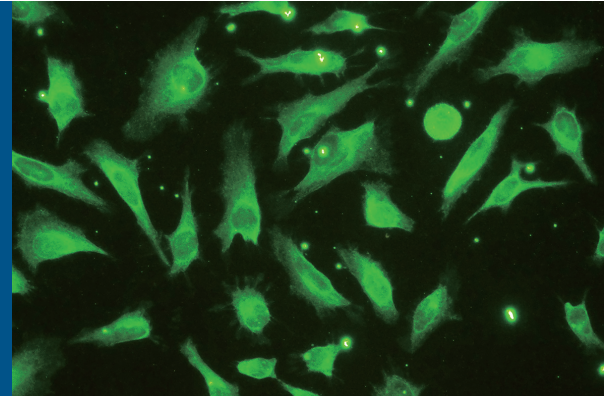

Features

LED light source

Wide-Field spectrum

Simple operation user interface

Low photo-bleaching rate

Reverse light output with liquid light guide

Custom designed thermal system for LED light engine

Benefits

Ecofriendly design, long life, instant on/off, stable light output, and no light bulb change.


Wide-Field Spectrum (350 nm ~ 700 nm), sufficient to use in all research occasions related to fluorescence, such as DAPI, GFP/FITC, mCherry/Texas Red, and Cy3 excitation.

Easy operation and programmable; clear status display on OLED display.

LED cold light used; increases the fluorescence cell sample survival rate.

Small light engine and liquid light guide design, perfect for research laboratories with limited space.

Custom designed thermal system effectively control the light engine/illuminator at the best operation temperature to ensure the illuminator performs at best conditions.

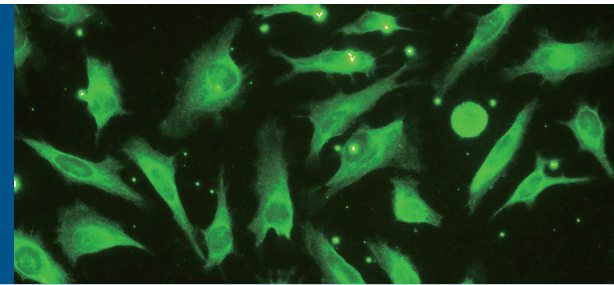


HYPER E301

Wide-Field Excitation Illuminator for Fluorescence Microscope

Specification

Wavelength range	350 ~ 700 nm
LED peaks/FWHM	365/20 nm, 460/40 nm, 560/80 nm
External power supply	Universal input 100-240 VAC, 50/60 Hz
Power consumption	180W
LED on/off response times	1 ms
Panel control	On/Off, 1% intensity increment
Dimensions (WxDxH) mm	128 x 338 x 210 (Excluding rubber legs, cap adaptor)
Weight	Approx. 5 Kg
Certifications	CE
Warranty (LEDs)	LEDs: 25,000 hours (or 3 years) Illuminator: 2 years (parts & labor)

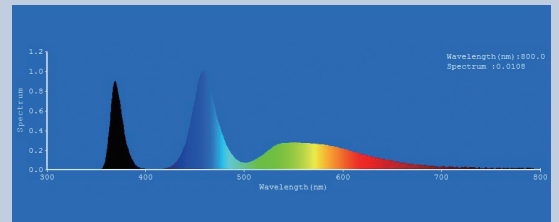


Accessories & Setup Illustration

Leica Nikon Olympus Zeiss

Collimator

Liquid Light Guide YODN Hyper E301 Illuminator

Dimensions

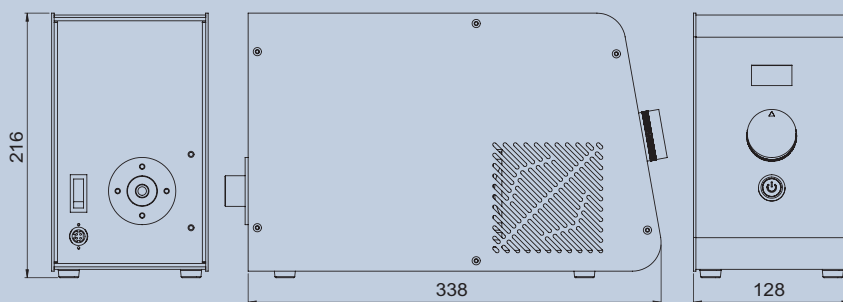

YODN Lighting Corp.

Address: 6F, No. 1, Creation Rd. II, Science-Based Industrial Park, Hsin-chu City, 30077 Taiwan

TEL: +886-3-5637218 FAX: +886-3-5794581

Contact: Ryan Shiau

Email: ryan_shiau@yodnlighting.com

What's APP ID: +886952968814



Follow Us on Facebook
www.facebook.com/YODN Lighting Corp.